



# EVO Series

HIGH SPEED FLOOR CENTRIFUGE

## EVO Series High Speed Floor Centrifuge

Versatile, high-throughput sample processing from 1 L bottles, to 50 mL conical tubes, to microtubes and PCR tubes per run. Save time and ensure sample integrity with fast acceleration and deceleration rates, combined with smooth runs at all speeds.

### Ergonomic Design

Touchscreen interface with bright, highly visible, and durable display simplifies run set-up.

### Adjustable Acceleration and Deceleration

15 rates available for to ensure the unit provide the best acceleration and deceleration profile to ensure your samples get the best separation.



### DISPLAY

Large Multi-color Touchscreen Display provide clear indication of running condition including rotor, speed, RCF, time etc.

The glove friendly touchscreen display gives you ease of control and a modern feel. An advanced microprocessor controls allow for setting/displaying of the rotational speed in either RPM or RCF. It is splash-proof and protects the inner electronics from external elements.

### ECO MODE



The display and control panel go into standby after a set time (adjustable), screen dims leaving only time on. Holds the chamber at a moderate temperature when idle to reduce compressor cycles while still enabling fast cool-down when restarted.

## PERFORMANCE

There are 15 rates available for customised acceleration and deceleration profile plus zero brake for each of your program. Ensuring your samples get the best separation. Frequency variable motor drive system enable the unit to operate in a stable and quiet manner. Pre-cool function allows set temperature to be reached rapidly.

## FLEXIBILITY

Our floor standing high speed centrifuges simplify high-throughput sample processing with up to **64,800 x g**, maximized capacity up to 6 liters, and an extensive rotor portfolio. Supporting labware from microtubes, to 1 Liter bottles. Overcome challenges of a shared laboratory setting with the productivity, safety, and reliability of a Newton Scientific Centrifuge

**Precision metal components for both shaft and rotor ensure reliable performance and long life**




## CONVENIENT

The EVO series high performance centrifuge enable the users to edit and store up to 99 programs with combination of speed, centrifugal force, running time, acceleration/ deceleration rate which is sufficient for all users in a lab.

Programmable Automatic lid release after run — The lid can be program to open once the rotor has come to a complete stop, saving time for the user.

Motorised lid lock\* - beneficial for users with limited physical strength or mobility, making it easier for anyone to operate the centrifuge without applying force when closing the lid.

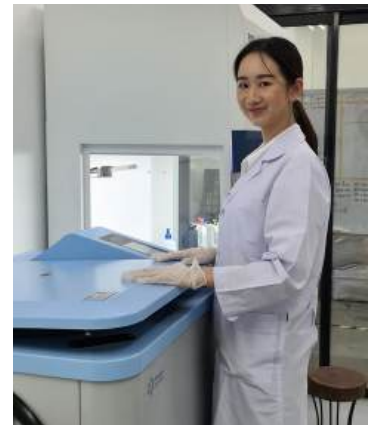

## SAFETY

Enhanced safety provided by electronic door lock and protection against imbalance, over speed and over heat. These three layers of protection ensure safety of the machine and user.

\*Available as standard on EVO 10 and EVO 25

## Ordering Information and Specifications



| Model          | EVO 23               | EVO 10               | EVO 25               |
|----------------|----------------------|----------------------|----------------------|
| Cat. No.       | 87002300             | 87001000             | 87002500             |
| Max RPM (rpm)  | 23,000               | 10,000               | 25,000               |
| RPM accuracy   | ±10 rpm              | ±10 rpm              | ±20 rpm              |
| Max RCF (x g)  | 35,500               | 18,300               | 64,800               |
| Max Capacity   | 4 x 750ml            | 6 x 1000ml           | 4 x 1000ml           |
| Set Time       | 1s to 23h 59m        | 1s to 99h 99m 99s    | 1s to 99h 99m 99s    |
| Temp Range     | -20°C to 40°C        | -20°C to 40°C        | -20°C to 40°C        |
| Temp Accuracy  | ±1.0°C               | ±1.0°C               | ±2.0°C               |
| Dimension (mm) | 660 x 520 x 820      | 905 x 745 x 1100     | 930 x 770 x 1080     |
| Weight (Kg)    | 140                  | 280                  | 300                  |
| Power Supply   | AC 200-240V, 50/60Hz | AC 200-240V, 50/60Hz | AC 200-240V, 50/60Hz |
| Power (W)      | 1300                 | 4500                 | 4500                 |
| Sound Level    | ≤ 65dBA              | ≤ 63dBA              | ≤ 60dBA              |



## Rotor and Accessories

EVO 23

| Cat. No. | Type    | Capacity   | Max. Speed (RPM) | RCF (x g) |
|----------|---------|--|------------------|-----------|
| 88200999 | Angle   | 12 x 1.5/2.0ml   | 23,000           | 38,371    |
| 88201000 | Angle   | 16 x 1.5/2.0ml   | 21,000           | 33,847    |
| 88201001 | Angle   | 24 x 1.5/2.0ml   | 18,000           | 31,127    |
| 88201003 | Angle   | 48 x 1.5/2.0ml   | 16,000           | 27,376    |
| 86005030 | Adaptor | 0.2 ml PCR tube (24 pcs)                                       |                  |           |
| 86005031 | Adaptor | 0.5 ml PCR tube (24 pcs)                                       |                  |           |
| 88201005 | Angle   | 12 x 5ml   | 18,500           | 23,199    |
| 88201006 | Angle   | 12 x 10ml  | 17,000           | 26,052    |
| 88213003 | Adaptor | 5 ml (12 pcs)  |                  |           |
| 88213004 | Adaptor | 1.5/2 ml (12 pcs)  |                  |           |
| 88201007 | Angle   | 12 x 15ml Con  | 13,000           | 19,185    |
| 88213005 | Adaptor | 10 ml (12 pcs)   |                  |           |
| 88213006 | Adaptor | 5 ml (12 pcs)  |                  |           |
| 88213007 | Adaptor | 1.5/2 ml (12 pcs)  |                  |           |
| 88201008 | Angle   | 6 x 50ml Con/Round<br>(includes 15ml Conical adaptors, 6 pcs ) | 14,000           | 22,106    |
| 88201009 | Angle   | 8 x 50ml Con/Round   | 12,000           | 14,148    |
| 88213008 | Adaptor | 15 ml Con (1pc)  |                  |           |
| 88213009 | Adaptor | 10 ml (1 pc)   |                  |           |
| 88213010 | Adaptor | 5 ml (1 pc)  |                  |           |
| 88213011 | Adaptor | 1.5/2 ml (1 pc)  |                  |           |
| 88201011 | Angle   | 6 x 100ml  | 12,500           | 19,085    |
| 88213012 | Adaptor | 50 ml (1 pc)   |                  |           |
| 88213013 | Adaptor | 15 ml (1 pc)   |                  |           |
| 88213014 | Adaptor | 2 x 10 ml (1 pc)   |                  |           |
| 88213015 | Adaptor | 3 x 5 ml (1 pc)  |                  |           |
| 88213016 | Adaptor | 4 x 1.5/2 ml (1 pc)  |                  |           |
| 87140250 | Angle   | 4 x 250ml  | 12,000           | 21,630    |
| 87140500 | Angle   | 4 x 500ml  | 8,000            | 10,375    |



**88201003**  
48 x 1.5/2ml



**88201008**  
6 x 50ml



**87140500**  
4 x 500ml

## Rotor and Accessories



### EVO 23

| Cat. No. | Type      | Capacity                            | Max. Speed (RPM) | RCF (x g) |
|----------|-----------|-------------------------------------|------------------|-----------|
| 88201017 | Swing Out | 4 x 500ml + Round Buckets           | 4,200            | 3,510     |
| 88213017 | Adapter   | 250 ml Flat Bottom (set of 4)       |                  |           |
| 88213018 | Adapter   | 2 x 100 ml Round Bottom (set of 4)  |                  |           |
| 88213019 | Adapter   | 4 x 50 ml Conical Bottom (set of 4) |                  |           |
| 88213020 | Adapter   | 4 x 50 ml Flat Bottom (set of 4)    |                  |           |
| 88213021 | Adapter   | 4 x 50 ml Round Bottom (set of 4)   |                  |           |
| 88213022 | Adapter   | 9 x 15 ml Conical Bottom (set of 4) |                  |           |
| 88213023 | Adapter   | 9 x 10 ml (set of 4)                |                  |           |
| 88213024 | Adapter   | 14 x 10 ml Blood Tube (set of 4)    |                  |           |
| 88213025 | Adapter   | 19 x 2-7 ml Blood Tube (set of 4)   |                  |           |
| 88213026 | Adapter   | 17 x 1.5/2 ml (set of 4)            |                  |           |
| 88201018 | Swing Out | 4 x 500ml + Round Buckets + Lids    | 4,200            | 3,510     |
| 88213017 | Adapter   | 250 ml Flat Bottom (set of 4)       |                  |           |
| 88213027 | Adapter   | 100 ml Round Bottom (set of 4)      |                  |           |
| 88213028 | Adapter   | 3 x 50 ml Conical Bottom (set of 4) |                  |           |
| 88213029 | Adapter   | 3 x 50 ml Flat Bottom (set of 4)    |                  |           |
| 88213030 | Adapter   | 3 x 50 ml Round Bottom (set of 4)   |                  |           |
| 88213022 | Adapter   | 9 x 15 ml Conical Bottom (set of 4) |                  |           |
| 88213023 | Adapter   | 9 x 10 ml (set of 4)                |                  |           |
| 88213031 | Adapter   | 13 x 10 ml Blood Tube (set of 4)    |                  |           |
| 88213032 | Adapter   | 17 x 2-7 ml Blood Tube (set of 4)   |                  |           |
| 88213026 | Adapter   | 17 x 1.5/2 ml (set of 4)            |                  |           |

## Rotor and Accessories



### EVO 23

| Cat. No. | Type      | Capacity                            | Max. Speed (RPM) | RCF (x g) |
|----------|-----------|-------------------------------------|------------------|-----------|
| 88201019 | Swing Out | 4 x 500ml + Rectangular Buckets     | 4,200            | 3,510     |
| 88213033 | Adapter   | 500 ml Flat Bottom (set of 4)       |                  |           |
| 88213034 | Adapter   | 250 ml Flat Bottom (set of 4)       |                  |           |
| 88213035 | Adapter   | 2 x 100 ml Round Bottom (set of 4)  |                  |           |
| 88213036 | Adapter   | 4 x 50 ml Conical Bottom (set of 4) |                  |           |
| 88213037 | Adapter   | 4 x 50 ml Flat Bottom (set of 4)    |                  |           |
| 88213038 | Adapter   | 4 x 50 ml Round Bottom (set of 4)   |                  |           |
| 88213039 | Adapter   | 9 x 15 ml Conical Bottom (set of 4) |                  |           |
| 88213040 | Adapter   | 9 x 15 ml Round Bottom (set of 4)   |                  |           |
| 88213041 | Adapter   | 16 x 10 ml Blood Tube (set of 4)    |                  |           |
| 88213042 | Adapter   | 29 x 2-7 ml Blood Tube (set of 4)   |                  |           |
| 88213043 | Adapter   | 17 x 1.5/2 ml (set of 4)            |                  |           |

## Rotor and Accessories



|          |           |  | EVO 23           |           |
|----------|-----------|--|------------------|-----------|
| Cat. No. | Type      | Capacity                                   | Max. Speed (RPM) | RCF (x g) |
| 88201020 | Swing Out | 4 x 750ml + Round Buckets                  | 4,200            | 3,777     |
| 88213044 | Adapter   | 500 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213045 | Adapter   | 500 ml Flat Bottom (set of 4)              |                  |           |
| 88213046 | Adapter   | 250 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213047 | Adapter   | 250 ml Flat Bottom (set of 4)              |                  |           |
| 88213048 | Adapter   | 4 x 100 ml Round Bottom (set of 4)         |                  |           |
| 88213049 | Adapter   | 7 x 50 ml Conical Bottom (set of 4)        |                  |           |
| 88213050 | Adapter   | 7 x 50 ml Flat Bottom (set of 4)           |                  |           |
| 88213051 | Adapter   | 7 x 50 ml Round Bottom (set of 4)          |                  |           |
| 88213052 | Adapter   | 14 x 15 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213053 | Adapter   | 14 x 10 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213052 | Adapter   | 14 x 15 ml Conical Bottom (set of 4)       |                  |           |
| 88213053 | Adapter   | 14 x 15 ml Round Bottom (set of 4)         |                  |           |
| 88213054 | Adapter   | 21 x 10 ml Blood Tube (set of 4)           |                  |           |
| 88213055 | Adapter   | 30 x 2-7 ml Blood Tube (set of 4)          |                  |           |
| 88213056 | Adapter   | 30 x 1.5/2 ml (set of 4)                   |                  |           |
| 88213057 | Adapter   | 250/300 ml Blood Bag (set of 4)            |                  |           |
| 88213058 | Carrier   | MTP/DWP (set of 4 pcs)                     |                  |           |

## Rotor and Accessories



|          |           |  | EVO 23           |           |
|----------|-----------|--|------------------|-----------|
| Cat. No. | Type      | Capacity                                   | Max. Speed (RPM) | RCF (x g) |
| 88201021 | Swing Out | Rotor, 4 x 750ml + Round Buckets + Lids    | 4,200            | 3,777     |
| 88213044 | Adapter   | 500 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213045 | Adapter   | 500 ml Flat Bottom (set of 4)              |                  |           |
| 88213046 | Adapter   | 250 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213047 | Adapter   | 250 ml Flat Bottom (set of 4)              |                  |           |
| 88213048 | Adapter   | 4 x 100 ml Round Bottom (set of 4)         |                  |           |
| 88213049 | Adapter   | 7 x 50 ml Conical Bottom (set of 4)        |                  |           |
| 88213050 | Adapter   | 7 x 50 ml Flat Bottom (set of 4)           |                  |           |
| 88213051 | Adapter   | 7 x 50 ml Round Bottom (set of 4)          |                  |           |
| 88213052 | Adapter   | 14 x 15 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213053 | Adapter   | 14 x 10 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213054 | Adapter   | 21 x 10 ml Blood Tube (set of 4)           |                  |           |
| 88213055 | Adapter   | 30 x 2-7 ml Blood Tube (set of 4)          |                  |           |
| 88213056 | Adapter   | 30 x 1.5/2 ml (set of 4)                   |                  |           |
| 88213057 | Adapter   | 450/500 ml Blood Bag (set of 4)            |                  |           |
| 88213058 | Carrier   | MTP/DWP (set of 4 pcs)                     |                  |           |

## Rotor and Accessories



### EVO 23

| Cat. No. | Type      | Capacity                                   | Max. Speed (RPM) | RCF (x g) |
|----------|-----------|--|------------------|-----------|
| 88201022 | Swing Out | 4 x 750ml + Rectangular Buckets            | 4,200            | 3,777     |
| 88213059 | Adapter   | 750 ml Flat Bottom (set of 4)              |                  |           |
| 88213060 | Adapter   | 500 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213061 | Adapter   | 500 ml Flat Bottom (set of 4)              |                  |           |
| 88213062 | Adapter   | 250 ml Corning Conical Tube (set of 4)     |                  |           |
| 88213063 | Adapter   | 250 ml Flat Bottom (set of 4)              |                  |           |
| 88213064 | Adapter   | 4 x 100 ml Round Bottom (set of 4)         |                  |           |
| 88213065 | Adapter   | 7 x 50 ml Conical Bottom (set of 4)        |                  |           |
| 88213066 | Adapter   | 7 x 50 ml Flat Bottom (set of 4)           |                  |           |
| 88213067 | Adapter   | 7 x 50 ml Round Bottom (set of 4)          |                  |           |
| 88213068 | Adapter   | 17 x 15 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213069 | Adapter   | 17 x 10 ml Conical/Round Bottom (set of 4) |                  |           |
| 88213070 | Adapter   | 28 x 10 ml Blood Tube (set of 4)           |                  |           |
| 88213071 | Adapter   | 34 x 2-7 ml Blood Tube (set of 4)          |                  |           |
| 88213072 | Adapter   | 30 x 1.5/2 ml (set of 4)                   |                  |           |
| 88213073 | Adapter   | 450/500 ml Blood Bag (set of 4)            |                  |           |
| 88213058 | Carrier   | MTP/DWP (set of 4 pcs)                     |                  |           |

## Cell Culture Packages



|          |   | EVO 23           |           |
|----------|---|------------------|-----------|
| Cat. No. | Capacity  | Max. Speed (RPM) | RCF (x g) |
| 88213075 | <b>Cell Culture Package 1:</b><br>1) Swing Out Rotor with rectangular buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 17 x 15 ml conical tube (set of 4). 68 tubes per run                                   | 4,200            | 3,777     |
| 88213076 | <b>Cell Culture Package 2:</b><br>1) AUTO LOCK Swing out rotor with rectangular buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 17 x 15 ml conical tube (set of 4). 68 tubes per run                         | 4,200            | 3,777     |
| 88213077 | <b>Cell Culture Package 3:</b><br>1) Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 14 x 15 ml conical tube (set of 4). 56 tubes per run   | 4,200            | 3,777     |
| 88213078 | <b>Cell Culture Package 4:</b><br>1) Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 14 x 15 ml conical tube (set of 4). 56 tubes per run<br>4) Sealing lids (set of 4)           | 4,200            | 3,777     |
| 88213079 | <b>Cell Culture Package 5:</b><br>1) AUTO LOCK Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 14 x 15 ml conical tube (set of 4). 56 tubes per run                               | 4,200            | 3,777     |
| 88213080 | <b>Cell Culture Package 6:</b><br>1) AUTO LOCK Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 7 x 50 ml conical tube (set of 4). 28 tubes per run<br>3) Adapter, 14 x 15 ml conical tube (set of 4). 56 tubes per run<br>4) Sealing lids (set of 4) | 4,200            | 3,777     |
| 88213081 | <b>Cell Culture Package 7:</b><br>1) Swing Out Rotor with buckets, 4 x 500ml<br>2) Adapter, 4 x 50 ml conical tube (set of 4). 16 tubes per run<br>3) Adapter, 9 x 15 ml conical tube (set of 4). 36 tubes per run  | 4,200            | 3,510     |
| 88213082 | <b>Cell Culture Package 8:</b><br>1) Swing out rotor with buckets, 4 x 500ml<br>2) Adapter, 3 x 50 ml conical tube (set of 4). 12 tubes per run<br>3) Adapter, 9 x 15 ml conical tube (set of 4). 36 tubes per run<br>4) Sealing lids (set of 4)                  | 4,200            | 3,510     |
| 88213083 | <b>Cell Culture Package 9:</b><br>1) Swing Out Rotor 4 x 750ml<br>2) Carrier for MTP/DWP  | 4,200            | 3,126     |

## Blood Tube Packages



|          |  | EVO 23           |           |
|----------|--|------------------|-----------|
| Cat. No. | Capacity   | Max. Speed (RPM) | RCF (x g) |
| 88213084 | <b>Blood Tube Package 1:</b><br>1) Swing Out Rotor with rectangular buckets, 4 x 750ml<br>2) Adapter, 34 x 2-7 ml blood tube (set of 4). 136 tubes per run<br>3) Adapter, 28 x 10 ml blood tube (set of 4). 112 tubes per run                                  | 4,200            | 3,777     |
| 88213085 | <b>Blood Tube Package 2:</b><br>1) AUTO LOCK Swing out rotor with rectangular buckets, 4 x 750ml<br>2) Adapter, 34 x 2-7 ml blood tube (set of 4). 136 tubes per run<br>3) Adapter, 28 x 10 ml blood tube (set of 4). 112 tubes per run                        | 4,200            | 3,777     |
| 88213086 | <b>Blood Tube Package 3:</b><br>1) Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 30 x 2-7 ml blood tube (set of 4), 120 tubes per run<br>3) Adapter, 21 x 10 ml blood tube (set of 4), 84 tubes per run   | 4,200            | 3,777     |
| 88213087 | <b>Blood Tube Package 4:</b><br>1) Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 30 x 2-7 ml blood tube (set of 4), 120 tubes per run<br>3) Adapter, 21 x 10 ml blood tube (set of 4), 84 tubes per run<br>4) Sealing lids (set of 4)           | 4,200            | 3,777     |
| 88213088 | <b>Blood Tube Package 5:</b><br>1) AUTO LOCK Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 30 x 2-7 ml blood tube (set of 4), 120 tubes per run<br>3) Adapter, 21 x 10 ml blood tube (set of 4), 84 tubes per run                               | 4,200            | 3,777     |
| 88213089 | <b>Blood Tube Package 6:</b><br>1) AUTO LOCK Swing Out Rotor with round buckets, 4 x 750ml<br>2) Adapter, 30 x 2-7 ml blood tube (set of 4), 120 tubes per run<br>3) Adapter, 21 x 10 ml blood tube (set of 4), 84 tubes per run<br>4) Sealing lids (set of 4) | 4,200            | 3,777     |
| 88213090 | <b>Blood Package 7:</b><br>1) Swing out rotor with buckets, 4 x 500ml<br>2) Adapter, 19 x 2-7 ml blood tube (set of 4), 76 tubes per run<br>3) Adapter, 14 x 10 ml blood tube (set of 4), 56 tubes per run   | 4,200            | 3,510     |
| 88213091 | <b>Blood Package 8:</b><br>1) AUTO LOCK Swing out rotor with buckets, 4 x 500ml<br>2) Adapter, 17 x 2-7 ml blood tube (set of 4), 68 tubes per run<br>3) Adapter, 13 x 10 ml blood tube (set of 4), 52 tubes per run<br>4) Sealing lids (set of 4)             | 4,200            | 3,510     |



|          |  | EVO 10           |           |
|----------|--|------------------|-----------|
| Cat. No. | Type                                   | Max. Speed (RPM) | RCF (x g) |
| 87160250 | NSA-6.250 Angle Rotor, 6 x 250ml       | 10,000           | 18,300    |
| 87160500 | NSA-6.500 Angle Rotor, 6 x 500ml       | 8,000            | 11,620    |
| 87161000 | NSA-6.1000 Angle Rotor, 6 x 1000ml     | 8,000            | 15,310    |
| 87241000 | NSS-4.1000 Swing-out Rotor, 4 x 1000ml | 4,200            | 4,430     |
| 87261000 | NSS-6.1000 Swing-out Rotor, 6 x 1000ml | 4,200            | 5,100     |

\*We provide custom solutions for specialise containers and vessels. Contact us with your requirement!



**87160500**  
6 x 500ml



**87161000**  
6 x 1000ml



**87261000**  
6 x 1000ml


**EVO 25**

| Cat. No. | Type   | Max. Speed (RPM) | RCF (x g) |
|----------|--|------------------|-----------|
| 87240002 | NSA-24.2T Angle Rotor, Titanium, 24×1.5/2.0ml  | 25,000           | 64,800    |
| 87180010 | NSA-18.10 Angle Rotor, 18×10ml   | 23,000           | 60,460    |
| 87180050 | NSA-8.50T Angle Rotor, 8×50ml Round Bottom, Titanium   | 21,000           | 52,356    |
| 87280050 | NSA-8.50 Angle Rotor, 8×50ml Round Bottom  | 20,000           | 46,420    |
| 87380050 | NSA-8.50C Angle Rotor, 8×50ml Conical  | 20,000           | 46,420    |
| 87480050 | NSA-8.50CF Angle Rotor, 8×50ml Conical, Carbon Fiber   | 20,000           | 46,420    |
| 87160050 | NSA-6.50 Angle Rotor, 6×50ml Conical/Round (includes 15ml Con adapters, 6 pcs)               | 18,000           | 33,506    |
| 87260050 | NSA-6.50 Angle Rotor, Carbon Fiber 6×50ml, Conical/Round (includes 15ml Con adapters, 6 pcs) | 18,000           | 33,506    |
| 87120015 | NSA-12.15C Angle Rotor, 12 x 15ml Conical  | 16,000           | 29,050    |
| 87160250 | NSA-6.250 Angle Rotor, 6 x 250ml   | 14,000           | 30,074    |
| 87260250 | NSA-6.250F Angle Rotor, Carbon Fiber, 6 x 250ml  | 14,000           | 30,074    |
| 87180100 | NSA-8.100 Angle Rotor, 8 x 100ml   | 13,000           | 19,830    |
| 87120050 | NSA-12.50 Angle Rotor, 6×50ml Con/Round  | 12,000           | 22,135    |
| 87160500 | NSA-6.500 Angle Rotor, 6 x 500ml   | 10,000           | 17,700    |
| 87260500 | NSA-6.500F Angle Rotor, Carbon Fiber 6 x 500ml   | 10,000           | 17,700    |
| 87141000 | NSA-4.1000 Angle Rotor, 4 x 1000ml   | 8,000            | 12,040    |
| 87241000 | NSA-4.1000F Angle Rotor, Carbon Fiber, 4 x 1000ml  | 8,000            | 12,040    |

\*We provide custom solutions for specialise containers and vessels. Contact us with your requirement!



**87141000**  
4 x 1000ml



**87160250**  
4 x 250ml



**87240002**  
24 x 1.5/2ml



Unit 2, New Horizon  
Business Centre, Barrows  
Road, Harlow, Essex, United  
Kingdom, CM19 5FN  
info@newton-scientific.com  
www.newton-scientific.com

